



Hydraulic & Engineering Instruments

Progressive in Design Perfect in Performance

OPERATING MANUAL



An ISO 9001 : 2015 Company



OPERATING INSTRUCTIONS

FOR

HS31.17 MOTORISED SIEVE SHAKER (WITH TIME SWITCH)

Hydraulic & Engineering Instruments
(An ISO 9001:2015 Company)

B-59/4, Naraina Industrial Area, Phase-II,

New Delhi-110 028 (INDIA)

Tel : 91-11-25893820 (4 Lines)

Fax : 91-11-25893152

E-mail : info@heicoin.com

Web : www.heicoin.com



ISO 9001:2015

C O N T E N T S

INTRODUCTION

DESCRIPTION

INSTALLATION

MAINTENANCE

INTRODUCTION

As a control on the raw materials or on the finished products dry sieve analysis are carried out regularly in a number of Industries. This can be used for ores, refractory materials, mineral aggregates, pigments, powdered coal, soap, cement, roofing materials, plastic molding powders, metal powders and in the pharmaceutical industry, in short wherever sieving is applied.

DESCRIPTION

The 'HEICO' Sieve Shaker is driven by a 1/2 H.P. Motor (3) through a reduction gear (4). The tapping base (9) is inclined from the vertical axis and the direction in addition to this gyratory motion of the tapping base (9) there is a stop pin (10) fixed under the tapping base (9). If the stop is removed the Shaker will have rotary motion also. The gear mechanism (5) runs in oil bath. A time switch (15) is also provided to facilitate setting the shaker to operate for required time. Suitable for operation in 220 V, 50 Cycles, single phase A.C. supply.

The gear oil (SAE 90 or equivalent) is supplied separately in a tin.

After unpacking the shaker, mount it on a firm level base. An oil level indicator (12), Oil inlet plug (11) and oil outlet plug (13) (painted Red) are provided. Fill the cylindrical drum Housing of the gear mechanism (5) with oil to the level indicated (12).

INSTALLATION


Nest the set of sieves together with lid and receiver and mount them on the tapping base.

Position the clamping plate (8) on top of the lid and tighten. Connect the supply lead to the main supply and set the point of the time switch (15) to the required time (Range of time switch is 0 to 60 minutes in 5 minutes divisions).

At the end of the interval the motor switches off automatically.

MAINTENANCE

Keep all the parts clean. Periodically remove the cover over the gear drive, connecting motor gear housing, and pack the pair of gears with high melting point grease.

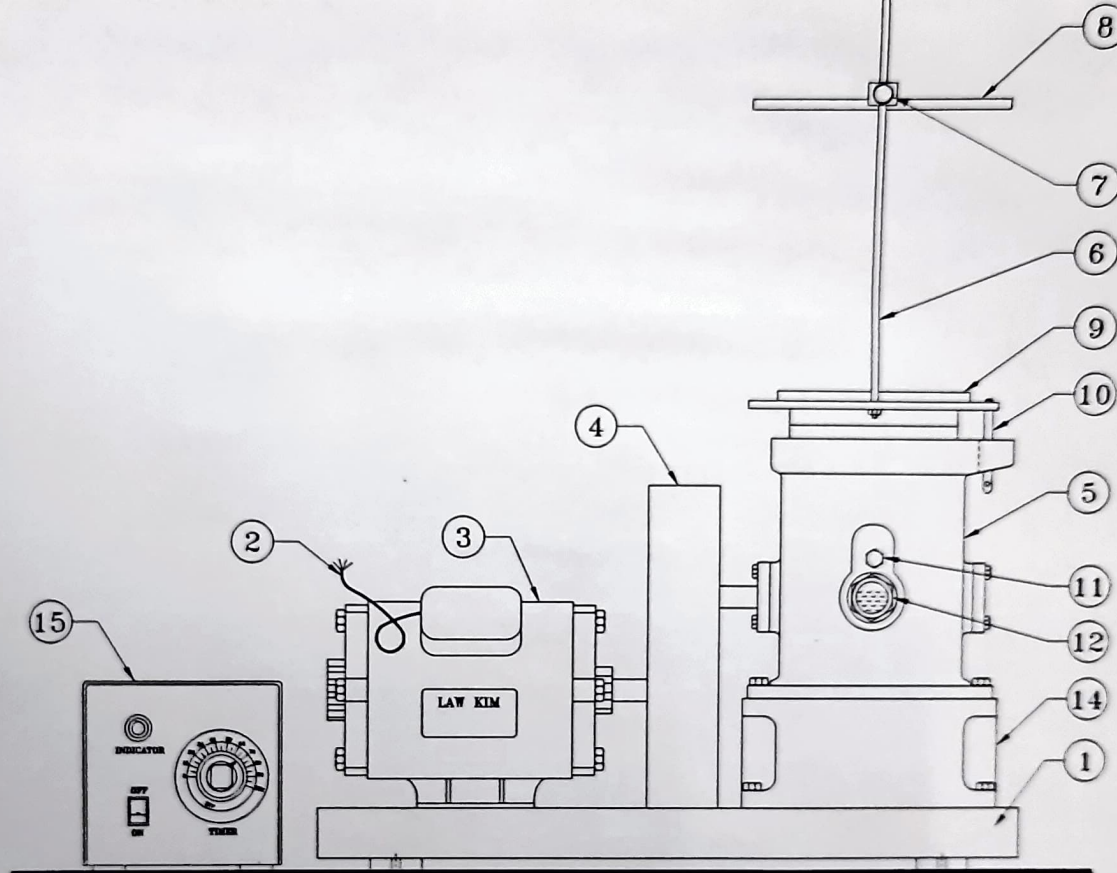

15	TIMER
14	SPACER
12	OIL LEVEL
11	OIL INLET/OUTLET
10	STOP PIN
9	TAPPING BASE
8	CLAMPING PLATE
7	LOCKING KNOB
6	SUPPORT RODS
5	GEAR MECHANISM
4	REDUCTION GEAR
3	MOTOR
2	MAIN LEAD
1	BASE
S.N. DESCRIPTION	

HEICO HYDRAULIC & ENGINEERING INSTRUMENTS A-13 NARAINA INDUSTRIAL AREA PHASE II NEW DELHI 110028 (INDIA) WEBSITE: www.heico.com PHONES 26060825-26		DRN KAMAL CHD H. SARDI APPD. SCALE: 1:1
MACHINE:- SIEVE SHAKER MOTORISED		CAT. NO. DRG. NO.
PART:- ASSEMBLY [HS31.17]		QTY. - REV. NO. -
MATERIAL:-		DR. - GR. -

THIS DRAWING IS THE PROPERTY OF HYDRAULIC AND ENGINEERING INSTRUMENTS AND SHOULD NOT BE REPRODUCED IN FULL OR PART THERE OF WITHOUT THEIR WRITTEN PERMISSION.